



MedicalDirector®

MedicalDirector Bluechip Day Surgery

December 2019 Update for South Australia



Health Communication Network Pty Limited Trading as MedicalDirector

ACN: 068 458 515 | ABN: 76 068 458 515 • Freecall: 1800 500 075 | Fax: +61 2 9460 7067

Level 5, 477 Pitt Street, Haymarket NSW 2000 • PO Box K1265 Haymarket, NSW 1240 Australia

We are pleased to announce the release of the December 2019 update for MedicalDirector Bluechip Day Surgery update for South Australia.

We strongly recommend that you apply the following update as this update is valid from December 2019.


If you have any issues, contact MedicalDirector Customer Support on 1300 300 161

Included in this Update

- Updated ICD10-AM tables to V11
- Updated HCP based extracts to reflect new ICD10 version
- Fixed problem with HCP extracts including morphology codes and/or leaving blanks in extract
- Updated National Private Patients Hospital Claim Form for 2019
- Updated Same Day Certificate for 2019 legislation
- Updated Grouper interface to support 3M CGS 2019.1 and DRG version 10
- Corrected errors in the ISAAC 1827 format extract
- Removed dropdown for year in ISAAC Control log
- Restored 4.0, 4.1, 4.2 as options for DRG Version in the Episodes form's Procedures tab
- Corrects error in the data export to ISAAC for episodes dated after 30/06/19. *Specifically, in field 146; "Contract provider hospital patient unit number"*


How to Backup Day Surgery Files

1. Download and extract the current SA Day Surgery file 2019-20 (DS_15.04_SA.zip) file to your local computer.
2. On your computer, locate and open the DSS folder. You can find this folder from within Bluechip, by selecting **Setup > Practice > Links > Day Surgery**
3. On your computer, within the DSS folder, create a 'backup' folder.
4. Now, go back up to the DSS folder. Locate and copy the old DSS, DSS Data, DSS Episodes and ISAAC Control Log files, and paste them into the **backup** folder.
5. Now, still in the **backup** folder, rename these files to indicate they're the 'old' e.g. DSS_old
6. Navigate back to where you downloaded and extracted the SA Day Surgery file 2019-20 (DS_15.04_SA.zip) file.
7. Copy the DSS.mdb file from the SA Day Surgery file 2019-20 (DS_15.04_SA.zip) folder and paste it into your DSS folder.

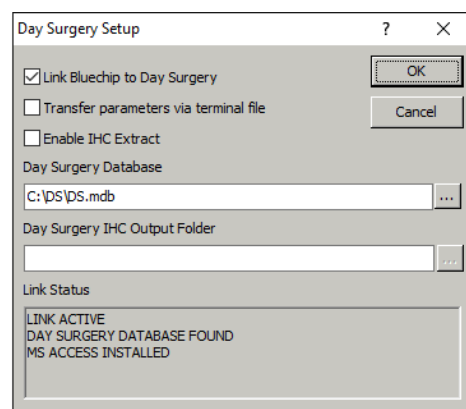

Relinking the Day Surgery module


1. Open Day Surgery. If you are prompted to create a new Episode, select **No**.

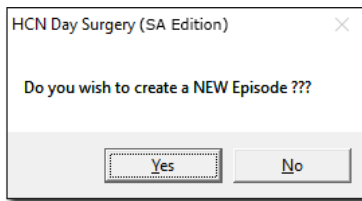

2. When prompted to **Find Episode**, select **Cancel**.



3. To re-link Episodes, select **Add-ins > Advanced > Relink Databases > Episodes**.
4. To re-link DS Data, select **Add-ins > Advanced > Relink Databases > DS Data**.
5. To confirm the re-linking you performed in steps 3 and 4, select **Add-ins > Advanced > Relink Databases > Show Links**

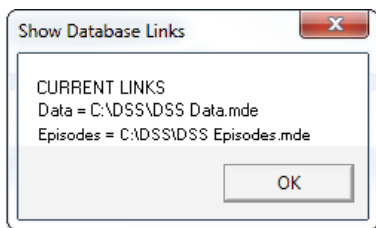

6. Select **Add-ins > Advanced > Settings** and update the DRG version to 9.0 or 10.0 if required.

Property	Setting	Comments
DefAttendingDoctor		Default Attending Doctor
defCountry	1100	Default value for Country
DisCancer	N	Disable cancer notification
DRGVersion	9.0	Default DRG Version

7. Close Day Surgery